

ASAHI FUBUKI Series

Neurovascular Guide Catheter

A high performance line of neurovascular

**ASAHI
FUBUKI**
Neurovascular Guide Catheter



Stable platform for your procedure.



**ASAHI
FUBUKI**
Neurovascular Guide Catheter Dilator Kit



A larger lumen guide with a smaller puncture site.

Neurovascular guide catheter docked with a dilator replacing the need for a sheath to minimize size of puncture site
*For this product line, a 'Stiff' type is available



Features of ASAHI FUBUKI 6,7,8 Fr and ASAHI FUBUKI Dilator Kit 4,5,6Fr

Flexibility

- Soft tip and flexible distal shaft minimizes vessel injury.

Pushability

- Well-balanced shaft with increased rigidity provides stable deployment and optimal pushability.

Lubricity

- Hydrophilic coating from the tip to 15cm (FUBUKI 8Fr, FUBUKI Dilator Kit 6Fr, and Stiff type = 5cm) allows smooth navigation through aortic arch and vessel ostium.

Visibility


- Radiopaque marker tip & shaft enhance visibility provided by a platinum marker on the tip and radiopaque material in the shaft.

FUBUKI Dilator Kit


*French sizes refer to the sheath size that would normally be required to introduce a guide catheter of the same size.

The 4/5/6 Fr FUBUKI Dilator Kit have the same shaft as the 6/7/8Fr regular FUBUKI, respectively.

ASAHI FUBUKI Dilator Kit (Sheath size)	I.D.	O.D.
4Fr	1.80mm (0.071")	2.09mm (6Fr)
5Fr	2.05mm (0.081")	2.40mm (7Fr)
6Fr	2.28mm (0.090")	2.70mm (8Fr)

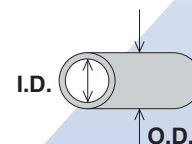

guide catheters in a wide range of sizes

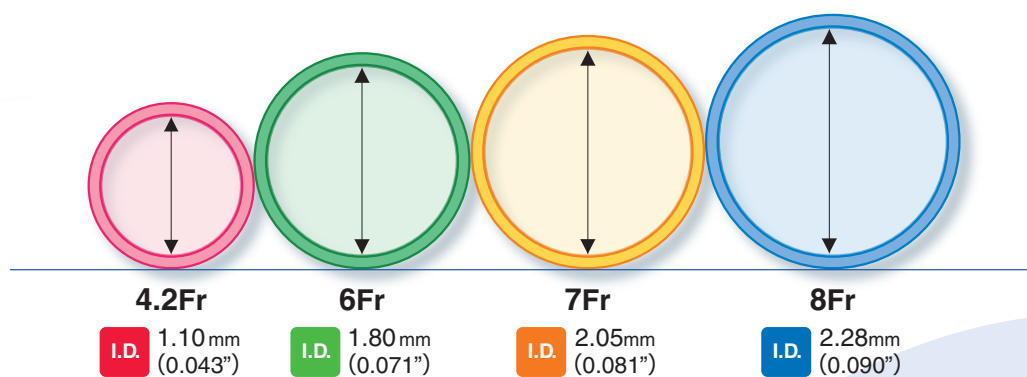


ASAHI
FUBUKI 043
Distal Support System




Strong support for delivering devices.

Features of ASAHI FUBUKI 043 4.2Fr ID: 1.10mm(0.043")

- Trackability**
 - Hydrophilic coating and flexible soft tube provides superb trackability in sharp bends or tortuous vessels.
- Enhanced kink resistance**
 - Proprietary wire braiding technology enhances kink resistance and maintains lumen integrity.
- Delivery Support**
 - Large inner lumen 1.10mm(0.043") with PTFE coating enables smooth delivery of the microcatheter.

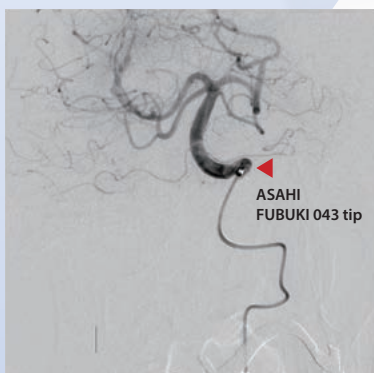


In meningioma treatment

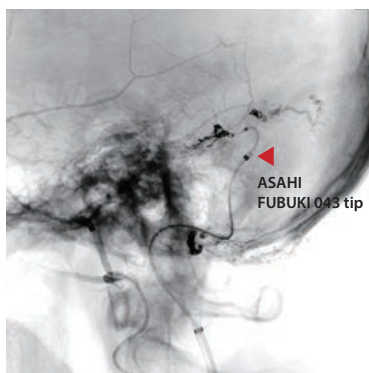

In AVF treatment

Structure & ordering information

■ Structure

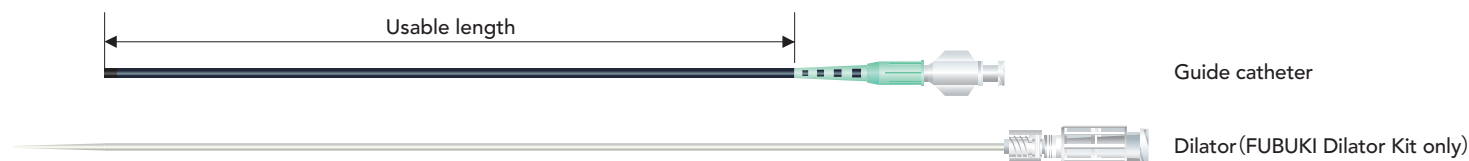

**ASAHI
FUBUKI**
Neurovascular Guide Catheter

Product	Catalog No.	Usable length	Coating length	Tip shape	Stiff type
ASAHI FUBUKI 6Fr	WAIN-FBK-6A	90 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-6AH	90 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-6S	90 _{cm}	15 _{cm}	Straight	—
	WAIN-FBK-6A80	80 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-6A80H	80 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-6S80	80 _{cm}	15 _{cm}	Straight	—
	WAIN-FBK-6AL	100 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-6ALH	100 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-6SL	100 _{cm}	15 _{cm}	Straight	—
	WAIN-FBK-6A110	110 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-6A110H	110 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-6S110	110 _{cm}	15 _{cm}	Straight	—
ASAHI FUBUKI 7Fr	WAIN-FBK-7A	90 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-7AH	90 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-7S	90 _{cm}	15 _{cm}	Straight	—
	WAIN-FBK-7A80	80 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-7A80H	80 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-7S80	80 _{cm}	15 _{cm}	Straight	—
	WAIN-FBK-7AL	100 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-7ALH	100 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-7SL	100 _{cm}	15 _{cm}	Straight	—
	WAIN-FBK-7A110	110 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-7A110H	110 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-7S110	110 _{cm}	15 _{cm}	Straight	—
ASAHI FUBUKI 8Fr	WAIN-FBK-8S	90 _{cm}	5 _{cm}	Straight	—
	WAIN-FBK-8S80	80 _{cm}	5 _{cm}	Straight	—
	WAIN-FBK-8SL	100 _{cm}	5 _{cm}	Straight	—
	WAIN-FBK-8S110	110 _{cm}	5 _{cm}	Straight	—

**ASAHI
FUBUKI**
Neurovascular Guide Catheter Dilator Kit

Product	Catalog No.	Usable length	Coating length	Tip shape	Stiff type
ASAHI FUBUKI Dilator Kit 4Fr	WAIN-FBK-4AD	90 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-4ADH	90 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-4SD	90 _{cm}	15 _{cm}	Straight	—
	WAIN-FBK-4AD80	80 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-4AD80H	80 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-4SD80	80 _{cm}	15 _{cm}	Straight	—
	WAIN-FBK-4ADL	100 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-4ADLH	100 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-4SDL	100 _{cm}	15 _{cm}	Straight	—
	WAIN-FBK-4AD110	110 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-4AD110H	110 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-4SD110	110 _{cm}	15 _{cm}	Straight	—
ASAHI FUBUKI Dilator Kit 5Fr	WAIN-FBK-5AD	90 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-5ADH	90 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-5SD	90 _{cm}	15 _{cm}	Straight	—
	WAIN-FBK-5AD80	80 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-5AD80H	80 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-5SD80	80 _{cm}	15 _{cm}	Straight	—
	WAIN-FBK-5ADL	100 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-5ADLH	100 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-5SDL	100 _{cm}	15 _{cm}	Straight	—
	WAIN-FBK-5AD110	110 _{cm}	15 _{cm}	Angled	—
	WAIN-FBK-5AD110H	110 _{cm}	5 _{cm}	Angled	●
	WAIN-FBK-5SD110	110 _{cm}	15 _{cm}	Straight	—
ASAHI FUBUKI Dilator Kit 6Fr	WAIN-FBK-6SD	90 _{cm}	5 _{cm}	Straight	—
	WAIN-FBK-6SD80	80 _{cm}	5 _{cm}	Straight	—
	WAIN-FBK-6SDL	100 _{cm}	5 _{cm}	Straight	—
	WAIN-FBK-6SD110	110 _{cm}	5 _{cm}	Straight	—

**ASAHI
FUBUKI 043**
Distal Support System

Product	Catalog No.	Usable length	Coating length	Tip shape
ASAHI FUBUKI 043	WAIN-FBK-4-120	120 _{cm}	105 _{cm}	Straight
	WAIN-FBK-4-125	125 _{cm}	110 _{cm}	Straight
	WAIN-FBK-4-130	130 _{cm}	115 _{cm}	Straight

Stylet and Peel-away included (FUBUKI 043 only)

CE 0344

© 2017 ASAHI INTECC CO., LTD.
"ASAHI" and "FUBUKI" are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries.

Manufactured by:

ASAHI INTECC

ASAHI INTECC CO.,LTD.

3-100 Akatsuki-cho, Seto, Aichi 489-0071 Japan

Contact phone number: +81-561-86-9101

<http://www.asahi-intecc.com>

Distributed by: