



Structure & ordering information

ASAHI Masters PARKWAY SOFT

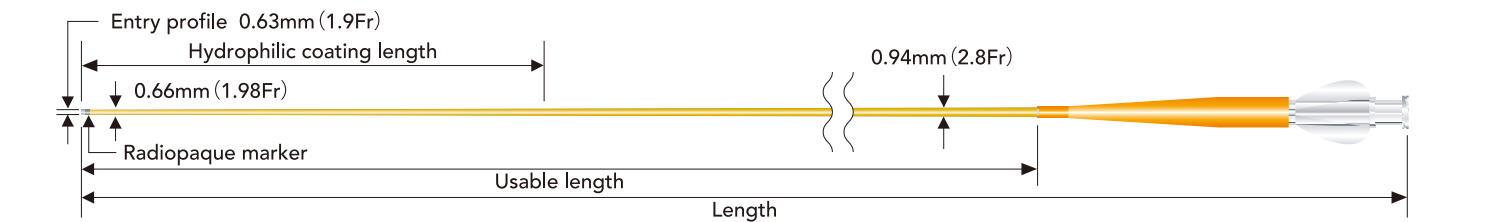


Product	Catalog No.	Usable length	Length	O.D. of distal shaft / O.D. of proximal shaft	Shaft I.D.	Max. guide wire	Hydrophilic coating length	Tip shape	Max. pressure
ASAHI Masters PARKWAY SOFT	WMST105-18PWSF	105cm	112cm	0.66mm (1.98Fr) 0.94mm (2.8Fr)	0.55mm (0.022inch)	0.46mm (0.018inch)	65cm	Straight	6895kPa (1000psi)
	WMST125-18PWSF	125cm	132cm	0.66mm (1.98Fr) 0.94mm (2.8Fr)	0.55mm (0.022inch)	0.46mm (0.018inch)	65cm	Straight	6895kPa (1000psi)
	WMST150-18PWSF	150cm	157cm	0.66mm (1.98Fr) 0.94mm (2.8Fr)	0.55mm (0.022inch)	0.46mm (0.018inch)	100cm	Straight	6895kPa (1000psi)
	WMST45A-18PWSSF	105cm	112cm	0.66mm (1.98Fr) 0.94mm (2.8Fr)	0.55mm (0.022inch)	0.46mm (0.018inch)	65cm	Angled 80°	6895kPa (1000psi)
	WMST45A-18PWSF	125cm	132cm	0.66mm (1.98Fr) 0.94mm (2.8Fr)	0.55mm (0.022inch)	0.46mm (0.018inch)	65cm	Angled 45°	6895kPa (1000psi)

ASAHI Masters PARKWAY HF

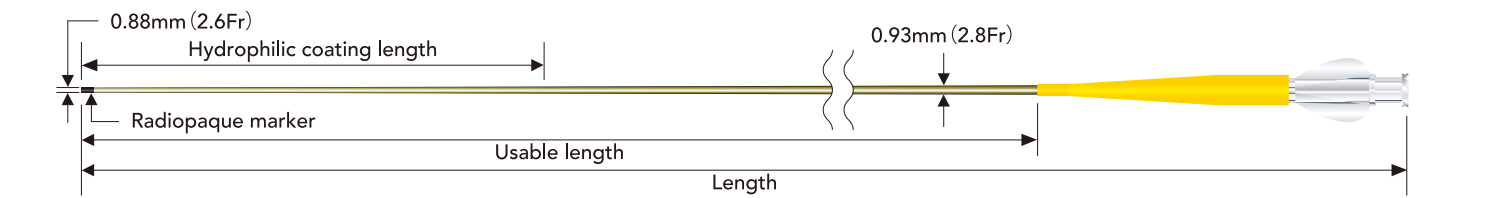

Product	Catalog No.	Usable length	Length	O.D. of distal shaft / O.D. of proximal shaft	Shaft I.D.	Max. guide wire	Hydrophilic coating length	Tip shape	Max. pressure
ASAHI Masters PARKWAY HF	WMST105-27HF	105cm	112cm	0.88mm (2.6Fr) 0.93mm (2.8Fr)	0.69mm (0.027inch)	0.53mm (0.021inch)	65cm	Straight	6895kPa (1000psi)
	WMST125-27HF	125cm	132cm	0.88mm (2.6Fr) 0.93mm (2.8Fr)	0.69mm (0.027inch)	0.53mm (0.021inch)	65cm	Straight	6895kPa (1000psi)
	WMST150-27HF	150cm	157cm	0.88mm (2.6Fr) 0.93mm (2.8Fr)	0.69mm (0.027inch)	0.53mm (0.021inch)	100cm	Straight	6895kPa (1000psi)

ASAHI Masters PARKWAY HF KIT

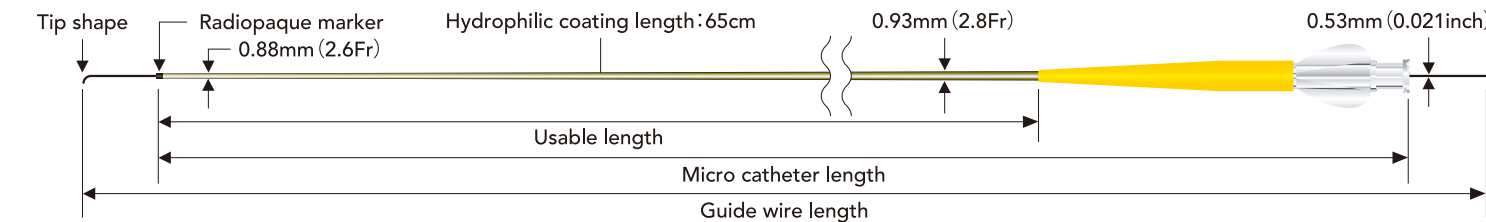

Product	Catalog No.	Micro catheter usable length	Micro catheter length	O.D. of distal shaft / O.D. of proximal shaft	Shaft I.D.	Guide wire shaft O.D.	Guide wire length	Tip shape of guidewire	Max. pressure
ASAHI Masters PARKWAY HF KIT	WMST105-27HFK	105cm	112cm	0.88mm (2.6Fr) 0.93mm (2.8Fr)	0.69mm (0.027inch)	0.53mm (0.021inch)	140cm	Angled 80°	6895kPa (1000psi)
	WMST125-27HFK	125cm	132cm	0.88mm (2.6Fr) 0.93mm (2.8Fr)	0.69mm (0.027inch)	0.53mm (0.021inch)	160cm	Angled 80°	6895kPa (1000psi)
	WMST105-27HFKM	105cm	112cm	0.88mm (2.6Fr) 0.93mm (2.8Fr)	0.69mm (0.027inch)	0.53mm (0.021inch)	140cm	Multi curve	6895kPa (1000psi)
	WMST125-27HFKM	125cm	132cm	0.88mm (2.6Fr) 0.93mm (2.8Fr)	0.69mm (0.027inch)	0.53mm (0.021inch)	160cm	Multi curve	6895kPa (1000psi)



© 2016 ASAHl INTECC CO., LTD.
"ASAHI" and "Masters PARKWAY" are trademarks or registered trademarks of ASAHl INTECC CO., LTD. in Japan and other countries.
Coated with SLIP-COAT® coating.

Manufactured by:



ASAHI INTECC CO.,LTD.

3-100 Akatsuki-cho, Seto, Aichi 489-0071 Japan

Contact phone number: +81-561-86-9101

<http://www.asahi-intecc.com>

Distributed by:

ASAHI Masters PARKWAY series/Ver.2/AMC-K16225

ASAHI Masters PARKWAY SOFT
ASAHI Masters PARKWAY HF
ASAHI Masters PARKWAY HF KIT


Microcatheter / Guide wire pre-loaded microcatheter



ASAHI Masters PARKWAY SOFT

Microcatheter

- Radius tip and small profile for distal accessibility.
- Large ID for high volume injection and smooth delivery of embolic materials.
- Tungsten braiding for lumen integrity and visibility.
- Pre-shaped tip.

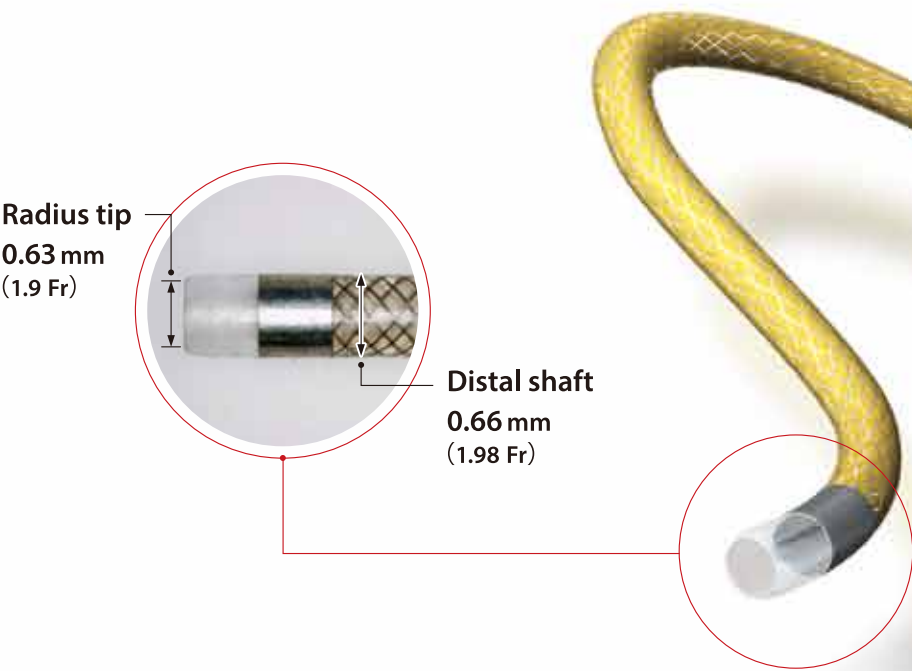

Flow rate

Product	Usable length	Contrast	Iodine content	Flow rate	Dead volume
ASAHI Masters PARKWAY SOFT	105 _{cm}	Iopamidol	300 (mg./ml)	2.6 (ml./sec)／6895 kPa (1000 psi)	0.43ml
			370 (mg./ml)	1.5 (ml./sec)／6895 kPa (1000 psi)	
	125 _{cm}		300 (mg./ml)	2.0 (ml./sec)／6895 kPa (1000 psi)	0.48ml
			370 (mg./ml)	1.0 (ml./sec)／6895 kPa (1000 psi)	
	150 _{cm}		300 (mg./ml)	1.6 (ml./sec)／6895 kPa (1000 psi)	0.54ml
			370 (mg./ml)	0.8 (ml./sec)／6895 kPa (1000 psi)	


ASAHI Masters PARKWAY HF

Microcatheter

ASAHI Masters PARKWAY HF KIT

Guide wire pre-loaded microcatheter

- Tungsten braiding for lumen integrity and visibility.
- High volume of contrast and therapeutic agents through 0.69 mm (0.027 inch) large inner lumen.
- 0.53 mm (0.021 inch) guide wire pre-loaded for Masters PARKWAY HF KIT.

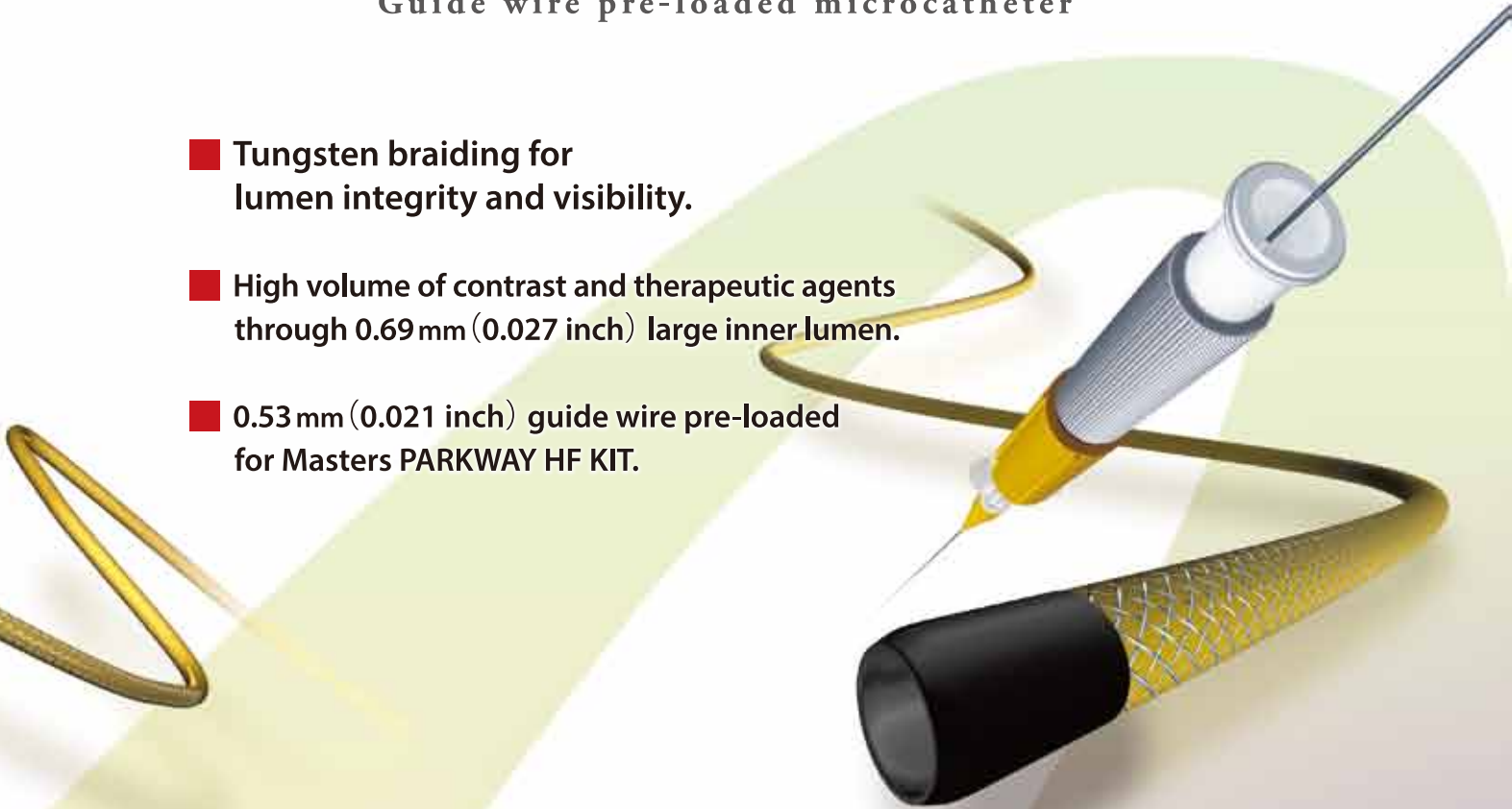

Flow rate

Product	Usable length	Contrast	Iodine content	Flow rate	Dead volume
ASAHI Masters PARKWAY HF	105 _{cm}	Iopamidol	300 (mg / ml)	5.1 (ml / sec) / 6895 kPa (1000 psi)	0.62 ml
			370 (mg / ml)	2.7 (ml / sec) / 6895 kPa (1000 psi)	
	125 _{cm}		300 (mg / ml)	4.4 (ml / sec) / 6895 kPa (1000 psi)	0.69 ml
			370 (mg / ml)	2.4 (ml / sec) / 6895 kPa (1000 psi)	
	150 _{cm}		300 (mg / ml)	3.8 (ml / sec) / 6895 kPa (1000 psi)	0.79 ml
			370 (mg / ml)	2.0 (ml / sec) / 6895 kPa (1000 psi)	